



PRODUCT-DATASHEET

SyriPIC[®]™ Needles

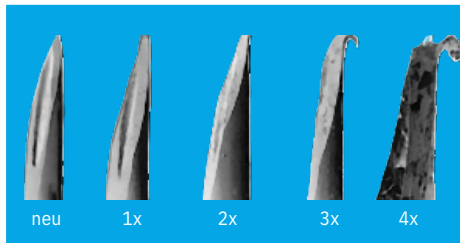


Safety

Many drug users face the challenge of having to get at their veins. A needle that glides through the skin as smoothly as possible and penetrates the vein safely is therefore of great benefit to these users in terms of safety and health.



The cutting technique is crucial to ensure that the needle can be used reliably.



A needle is generally intended for single use and should then be disposed of properly.

Different cutting techniques significantly affect the quality of a needle. A low-quality needle is damaged by repeated puncture attempts. These can blunt it and cause it to become unusable or very dangerous for the user.

SyriPIC[®]™ needles have a particularly sharp lancet tip. The exceptionally precisely-cut edges penetrate easily into the tissue and ensure a low-pain puncture. Silicone treatment of the tip allows the needle to penetrate gently.

**SyriPIC[®]™
Needles**









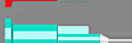
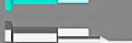
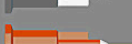




The specially thin-walled SyriPIC[®]™ needles allow for an increased rate of flow of the injection solution.


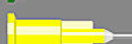

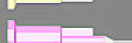
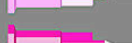

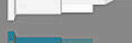

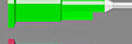
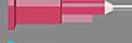




This means that needles with a reduced outer diameter can be used.

The high-quality materials and workmanship of the SyriPIC[®]™ needles are worth protecting effectively until they are used. This is achieved by means of a sturdy protective cap on the needle, which is made of resistant polypropylene. This is also transparent in order to make the ISO 6009 color-coded luer hub visible and thus ensure increased user safety.

Color Coding (ISO 6009)

	OD* (mm)	ID (mm)	Gauge	Colour**
	0.18	0.08	34G	Orange
	0.20	0.11	33G	Black
	0.23	0.11	32G	Deep Green
	0.25	0.12	31G	White
	0.30	0.16	30G	Yellow
	0.33	0.18	29G	Red
	0.36	0.18	28G	Blue-green
	0.40	0.21	27G	Medium Grey
	0.45	0.25	26G	Brown
	0.50	0.26	25G	Orange
	0.55	0.30	24G	Med Purple
	0.60	0.33	23G	Deep Blue
	0.70	0.41	22G	Black

	OD* (mm)	ID (mm)	Gauge	Colour**
	0.80	0.52	21G	Deep Green
	0.90	0.60	20G	Yellow
	1.10	0.67	19G	Cream
	1.20	0.84	18G	Pink
	1.40	1.07	17G	Red-violet
	1.60	1.25	16G	White
	1.80	1.40	15G	Blue-grey
	2.10	1.60	14G	Pale Green
	2.40	1.80	13G	Purple
	2.70	2.15	12G	Pale Blue
	3.00	2.40	11G	Green-yellow
	3.40	2.70	10G	Olive-brown

*ISO 9626 (Outer Diameter, OD)
**ISO 6009 (Colour)



**Injection needles
of the highest quality**
with particularly sharp lancet tips
and safety needles

SyriPIC[®]™ Needles

INFORMATION

Product Description

- All SyriPIC[®]™ needles from Compet Medical AG
 - comply with the latest injection technology
 - have a special thin wall
 - are made of stainless chrome-nickel steel
 - are subject to systematic checks
 - are individually sterile-packaged in a blister pack
 - operate at the highest level of quality.
- Low-pain & low-injury applications are ensured by
 - special bevelled long cut
 - puncture needles made of triple-tapered, stainless steel
 - extremely smooth metal surfaces - with silicone coating
 - lancet tips designed for exceptional sharpness.
- Transparent needle protection cap
 - with optimal needle protection prior to use
 - made of polypropylene.



Did you know that ...

... our SyriPIC[®]™ needles have another, very special advantage?

The specially thin wall.

*This makes it possible to choose a smaller needle size,
as the flow rate of the injection solution increases.*

Storage & Shipping Informationen

- Inner packaging:
 - Boxes of 100 pieces
 - Dimensions - box (L x W x H): 9.5 x 8.5 x 9.0 cm
- Sales unit:
 - 50 boxes of 100 pieces = 1 sales unit
 - up to max. 3 sales units per consignment



Smart ordering & cost-effective shopping:
If you bundle several sales units into one order,
shipping costs are reduced and you save money!

Product Portfolio

Art.-No. REF	Article Description	Outside Diameter		Length		Colour Code (ISO 6009)	Pieces per Sales Unit	Weight per Sales Unit (kg)	Max. Number of Sales Units per Shipment
		Gauge	mm	Inch	mm				
20520	SyriPIC [®] ™ - Needle 23G x 1"	23	0.60	1	25	blau	5'000	7.2	3
20420	SyriPIC [®] ™ - Needle 24G x 1"	24	0.55	1	25	violett	5'000	7.2	3
20220	SyriPIC [®] ™ - Needle 25G x 5/8"	25	0.50	5/8	16	orange	5'000	7.2	3
20120	SyriPIC [®] ™ - Needle 26G x 1/2"	26	0.45	1/2	12	braun	5'000	7.2	3
20320	SyriPIC [®] ™ - Needle 26G x 1"	26	0.45	1	25	braun	5'000	7.2	3
20720	SyriPIC [®] ™ - Needle 27G x 3/4"	27	0.40	3/4	19	grau	5'000	7.2	3
20820	SyriPIC [®] ™ - Needle 30G x 1/2"	30	0.30	1/2	12	gelb	5'000	7.2	3

Product Notes



SyriPIC[®]™ needles are medical devices.



The CE marking guarantees the consumer (in the EEA) uniform protection in health, safety and environmental matters.



Follow instructions for use.



Do not use products with damaged packaging.



Items are packaged sterile.



Products are intended for single use only.



Production lot, production & expiration date are marked on the packaging.



Storage: dry, cool and protected from the sun.

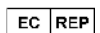
Production, Sales, Monitoring




Ningbo Greetmed Medical Instruments Co., Ltd.
16F-a, Building 1, No. 98 Chuangyuan Road, Hi-Tech Zone
315042 Ningbo, Zhejiang Prov. China



COMPET MEDICAL AG
Höhenweg 2b, 8834 Schindellegi, Switzerland



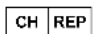Shanghai International Holding Corp. GmbH
Eiffelstrasse 80, 20537 Hamburg, Germany



Qunique GmbH
Bahnhofweg 17, 5610 Wohlen, Switzerland



Where there is hope, there is a future.

This fundamental conviction is what drives us at COMPET MEDICAL AG. Thus, we provide a wide range of user-friendly and high-quality products for drug assistance and addiction prevention (Harm Reduction).

The aim is to reduce health damage for drug users and to make the application as easy and safe as possible - to offer the best possible assistance for survival.

COMPET MEDICAL AG
Headquarters Switzerland
Höhenweg 2b
CH-8834 Schindellegi
+41 71 680 00 60

COMPET MEDICAL AG
Warehouse EU
Max Stromeyer Str.170
DE-78467 Konstanz
+49 7531 361 07 90

info@competmedical.com
www.competmedical.com



**COMPET
MEDICAL AG**

PRODUCT- DATASHEETS